

Multi-Function Examination Tables and Accessories



Distributore Ufficiale per l'Italia

REF ITALIA S.R.L

Via Nomentana, 403 – 00162 ROMA

Tel. +39 0645547880 – Fax +39 0645547881

Mail: info@refitalia.it

www.refitalia.it – www.colposcopi.it



Behandlung



Multi-Function Examination Tables

Multi-Function Examination Tables AGA-UNI-MED and AGA-DUO-LIFTO-MED

With their divided tabletop, comprehensive range of equipment options and included accessories, the multi-function examination tables can be used not only in general medicine but also in surgical, gynecological, proctology and rectoscopy practices.

The multi-function examination tables AGA-UNI-MED (Article No. A-2065/EE) and AGA-DUO-LIFTO-MED (Article No. H-U 2065/3EE) impress by their functionality, characteristics and durability.



1 A-2065/EE



2 H-U 2065/3EE

	AGA-UNI-MED	AGA-DUO-LIFTO-MED
Article number	A-2065/EE	H-U 2065/3EE
Picture(s)	1	2
Price		
Frame	Chrome	
Width (mm)	650 (710 with standard rails)	
Length (mm)	1900	
Head/Main/Foot section	500/1050/350	
Height adjustment (mm)	650-850 (Electric)	600-850 (Electric)
Trendelenburg/ Reverse Trendelenburg	25° (Electric)	15° (Electric)
Adjustment mechanism	Foot switches for simultaneous height and tilt adjustment	Hand module with 6 buttons for simultaneous height and tilt adjustment
Tabletop	3 sections	
Upholstery thickness (mm)	64	
Maximum load (kg)	200	150
Headrest adjustment	+40° (Clamping rods)	
Seat adjustment	+30° (Clamping rods)	
Backrest adjustment	+30° (Gas pressure spring)	
Pelvic cutout	✓	
Levelling feet	✓	
Paper roll holder	✓	
Accessories	Pages 14, 90-91	

Additional Features

- AGA-UNI-MED**
 The electric height and tilt adjustments are done with maintenance-free, spill-proof and extremely quiet electric motors, controlled by air pulses via foot switches (or a hand module, at no extra charge). Drive specifications (each motor): 230V, 50/60Hz, 3,7A, 850W, IPX4, v=14mm/s, 7000N.
- AGA-DUO-LIFTO-MED**
 The quiet lifting columns (low-voltage electric drives), made of anodized light alloy, are controlled via a transformer by a hand module (Multitasking). Drive specifications (each motor unit): 230V, 50/60Hz, 6,0A, 240W, IPX6, v=16mm/s, 2x1000N.

Both multi-function examination tables include:

- A pair of lateral standard rails (1720 mm),
- A pair of leg rests with plastic shells,
- A pair of mounting clamps,
- A stainless steel pan (300x200x65 mm),
- A paper roll holder,
- The foot switches/hand module.

Options

- The AGA-UNI-MED and AGA-DUO-LIFTO-MED are also available at extra charge with a black electrically conductive synthetic leather upholstery and a potential equalization socket (see page 4).

Multi-Function Examination Tables

Combined Multi-Function Tables with Height Adjustment

To prevent rocking and swinging movements of the tabletop, these combined multi-function tables are manufactured with precision dual-box rectangular profiles. They are available with manual, hydraulic oder electric height adjustment and offer a wide range of applications. The height adjustment take place without lateral offset and with little effort, even with heavier patients (up to 200 kg).



1 HK-2065/3H



2 HK-2065/3E

Article Number	HK-2065/3H	HK-2065/3Y	HK-2065/3E
Picture(s)	1	Similar 1, 2	2
Price			
Frame (tabletop/leg)	Chrome		
Width (mm)	650		
Length (mm) Head/Main/Foot section	1900 500/1050/350		
Height adjustment (mm)	650-850 (Manual)	650-850 (Hydraulic)	650-850 (Electric)
Tabletop	3 sections		
Upholstery thickness (mm)	64		
Maximum load (kg)	200		
Headrest adjustment	+40° (Clamping rods)		
Seat/backrest adjustment	+30° (Clamping rods)		
Pelvic cutout	✓		
Levelling feet	✓		
Rectoscopy bolster (Pic. 3)	RRH-65 *Retrofitting is not possible		
Paper roll holder	P-2065/A		
Accessories	Pages 90-91		

The combined multi-function tables with height adjustment include:

- A pair of leg rests with plastic shells,
- A pair of mounting clamps,
- A stainless steel pan (300x200x65 mm).

Options

- In order to treat patients in the knee-elbow position, the combined multi-function tables are also available with a four-piece tabletop and integrated, continuously adjustable rectoscopy bolster (Article No. RRH-65).



3 RRH-65

Additional Features

- **HK-2065/3H**
The manual height adjustment is done using maintenance-free gear with optimal transmission ratio (i=10mm/rev) for smooth, quick and easy operation of the hand crank from the left or the right-hand side.
- **HK-2065/3Y**
The hydraulic height adjustment is done with a centrally located

foot pump. The foot pedal can be mounted on either side (please specify upon ordering). A second foot pedal is available at extra charge.

- **HK-2065/3E**
The electric adjustment is done with a maintenance-free, extremely quiet and spill-proof electric motor, controlled by air pulses via a hand switch (or foot switch). Drive specifications: 230V, 50/60Hz, 3,7A, 850W, IPX4, v=14mm/s, 7000N.

Multi-Function Examination Tables

Combined Multi-Function Tables and AGA-MULTIFLEX Tables

The combined multi-function tables are also available with a fixed height, as presented hereafter. The seat upholstery and frame feature a pelvic cutout. These multi-function tables can therefore be used in the field of general medicine as well as for gynaecology and proctology (lithotomy position) examinations.

The AGA-MULTIFLEX tables were developed as an alternative to the combined multi-function tables. On these models, the folding leg section enables variable table length (from 1920 to 1380 mm) for more versatility of use.



1 K-1065/3



2 MULTIFLEX II



3 MULTIFLEX I



4 RRH-65

Options

- In order to treat patients in the knee-elbow position, the combined multi-function tables are also available with a four-piece tabletop and integrated, continuously adjustable rectoscopy bolster (Article No. RRH-65).

	Combined Multi-Function Tables			AGA-MULTIFLEX	
	K-1065/3	K-2065/3	K-3065/3	MULTIFLEX I	MULTIFLEX II
Article number	K-1065/3	K-2065/3	K-3065/3	MULTIFLEX I	MULTIFLEX II
Picture(s)	1	Similar 1	Similar 1	3	2
Price					
Frame	Lacquered	Chrome	Stainl. steel	Lacquered	
Length (mm)	1900			1900	
Head/Main/Foot section	500/1050/350			580/720/620	
Width (mm)	650			650	
Height (mm)	650/700/750/800			800	
Tabletop	3 sections			3 sections	
Upholstery thickness (mm)	64			64	
Maximum load (kg)	200			200	
Headrest adjustment	+40° (Clamping rods)			+30° (Grid)	
Seat/backrest adjustment	+30° (Clamping rods)			—	
Foot section adjustment	—			-90°	
Pelvic cutout	✓			✓	
Rectoscopy bolster	RRH-65 *Retrofitting not possible			—	
Leg rests (set of two)	✓			—	✓
Fixing clamps (set of two)	✓			—	✓
Stainless steel pan	✓			—	✓
Paper roll holder	P-2065			✓	
Accessories	Pages 90-91			Pages 90-91	

✓ included in delivery — not available

Multi-Function Examination Tables

Combined Multi-Function Tables AGA-BOZ

With its multiple configuration options, the combined multi-function examination table AGA-BOZ is both space-saving and versatile. Its main feature is the variable tabletop length, which can be shortened from 2100 mm to 1400 mm in a few steps. Due to the mechanism for sliding it under the tabletop for length adjustment, the leg section (width 530 mm) is slightly narrower than the main surface and the headrest (width 700 mm).

The two white melamine cabinets (594x594x318 mm, WxDxH) on models BOZ-1070/3S and BOZ-1070/3NS may be positioned individually as desired with the front panels facing forward/backward or to the right/left-hand side.



1 BOZ-1070/3S




2 BOZ-1070/3NS

Article Number	BOZ-1070/3S	BOZ-1070/3NS	BOZ-1070/3	BOZ-1070/3N
Picture(s)	1	2	3, 4	5
Price				
Frame	Lacquered			
Length (mm) Head/Main/Foot section	2100 500/800/800			
Width (mm)	700			
Height (mm)	800			
Tabletop	3 sections			
Upholstery thickness (mm)	64			
Maximum load (kg)	200			
Headrest adjustment	+30° (Grid)			
Paper roll holder (rod)	✓			
Two drawer units*	✓	✓	–	–
Step stool*	✓	–	✓	–
Pivoting stainless steel pan*	–	–	✓	–
Leg rests with plastic shells and fixing clamps (set of two)*	✓	–	✓	–
Accessories	Pages 90-91			

* Retrofitting is not possible

✓ included in delivery – not available

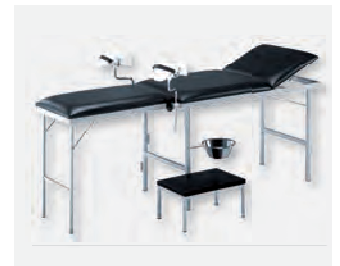

3 BOZ-1070/3



4 BOZ-1070/3













5 BOZ-1070/3N

Accessories

Accessories for Multi-Function Examination Tables

Picture	Article Number	Description	Price
	RLM-20	Knee roll, Ø 180x560 mm, synthetic leather, with zipper closure. Upholstery colours on page 4.	
	RM-12	Neck roll, Ø 140x400 mm, synthetic leather, with zipper closure. Upholstery colours on page 4.	
	KKM-16	Wedge-shaped headrest, 450x400x80/10 mm (LxWxH), synthetic leather, with zipper closure. Upholstery colours on page 4.	
	HRM-14	Half-round neck bolster, Ø 200x400 mm, synthetic leather, with zipper closure. Upholstery colours on page 4.	
	HRM-14/Band HRM-14/Klett	Half-round neck bolster with strap. Half-round neck bolster with hook-and-loop fasteners. Synthetic leather, Ø 200x400 mm. Upholstery colours on page 4.	
	P-2065 P-2080 P-20100	Paper roll holder, chrome-plated, width 650 mm. Paper roll holder, chrome-plated, width 800 mm. Paper roll holder, chrome-plated, width 1000 mm. Suitable for mounting on the leg frame of AGA tables.	
	P-2065/A P-2080/A	Paper roll holder, chrome-plated, with mounting hardware, width 650 mm. Paper roll holder, chrome-plated, with mounting hardware, width 800 mm. Suitable for mounting on the tabletop frame of the height-adjustable AGA multi-function tables.	
	LR-75 LR-125	Set of castors, Ø 75 mm. Set of castors, Ø 125 mm. A set consists of four castors, two of which with locking function. At time of ordering, please indicate the required total height of the table (including castors).	
	LR/F-75 LR/F-125	Mobile base on retractable castors, Ø 75 mm. Mobile base on retractable castors, Ø 125 mm. The foot pedal to secure/release is located at the head and foot end (press twice). At time of ordering, please indicate the required total height of the table (including castors). *Retrofitting is not possible.	
	LR/Z-125	Central locking castor system, Ø 125 mm, offering excellent mobility and comfort of use. All four double castors are connected together by rods. A foot pedal to secure/release is located at the head and foot ends of the table (press once). At time of ordering, please indicate the required total height of the table (including castors). *Retrofitting is not possible.	
	KVG	Headrest adjustment by gas pressure spring instead of grid or clamping rod mechanism, for most AGA multi-function tables. *Retrofitting is not possible.	

Accessories for Multi-Function Examination Tables

Picture	Article Number	Description	Price
	GS-1710 EfGS	Face cutout in the headrest, oval shape, 205x105 mm (LxW). *Retrofitting is not possible. Upholstered insert to close the face cutout.	
	PG-220 PG-260	Patient strapping belt, dimensions 2200x250 mm, LxW. Patient strapping belt, dimensions 2600x250 mm, LxW. To secure the patient's position. Upholstery colours on page 4.	
	OPS-40/A OPS-60/A OPS-100/A OPS-140/A OPS-170/A	Standard rails (25x10 mm, LxW), set of two, stainless steel, with spacers. Available in different lengths (the digits in the article number indicate the respective length in cm). Please specify the desired position of the rail on the table frame (head, foot or center of the table) at time of order. *Retrofitting is possible to a limited extent only.	
	SG-1000 SG-2000 SG-3000	Side guard, lacquered, dimensions 800x250 mm (LxH). Side guard, chrome-plated, dimensions 800x250 mm (LxH). Side guard, stainless steel, dimensions 800x250 mm (LxH). *Retrofitting of the removable side guards is not possible.	
	INF-1720/VA INF-30170/BKO	Infusion (IV) holder stainless steel, pole length 1700 mm, with four hooks and with fixing clamp (screw system). Please specify the exact desired position of the fixing clamp on the tabletop frame. *Retrofitting is possible to a limited extent only. Infusion (IV) holder, stainless steel, length 1700 mm, with four hooks and a mounting clamp. *Standard rail(s) required.	
	BKG-600	Göpel mounting clamps (set of two), adjustable, chrome-plated, suitable for standard rails. *Standard rails required.	
	BHG-600	Göpel leg rests (set of two), with spherical joints, PU-foam pads with black upholstery. *Mounting clamps (set of two) for standard rails required.	
	RRH-65	Rectoscopy bolster, continuously height-adjustable, dimensions 740x65x550 mm (LxWxH). *Mounting clamps (set of two) required.	
	RL-2060/BK	Ring armrests (set of two), black, dimensions 270x340x45 mm (HxLxW). *Mounting clamps for standard rails included.	
	A-2047	Arm support, chrome-plated, dimensions 400x120 mm (LxW), with swivel mechanism, PU-foam pad with black upholstery. *Mounting clamp for standard rail included.	