

## Universal Tables and Accessories

---



Distributore Ufficiale per l'Italia

**REF ITALIA S.R.L**

Via Nomentana, 403 – 00162 ROMA

Tel. +39 0645547880 – Fax +39 0645547881

Mail: [info@refitalia.it](mailto:info@refitalia.it)

[www.refitalia.it](http://www.refitalia.it) – [www.colposcopi.it](http://www.colposcopi.it)





# Universal Tables

## Examination Units

The examination units presented here, consisting of a universal table with (mobile) drawer units and optionally a swivel stool, are available as sets at a competitive price.


These products are practical combinations for a quick and easy setup of examination rooms.



1 BUD-6265



2 SSL-6558

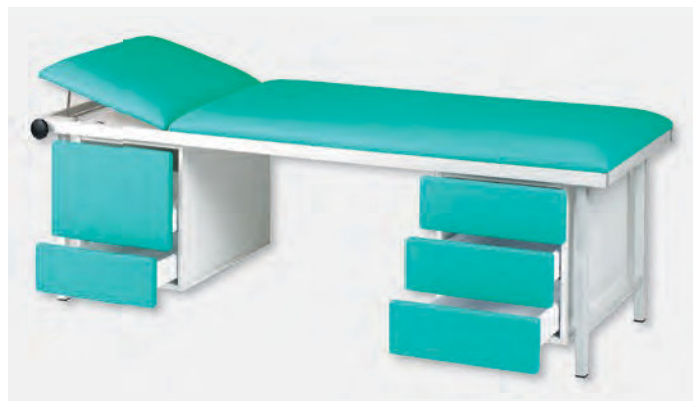

3 UEH-5865

### BUD-6265

This unit includes an examination table (similar to ULI-6000, page 45) with paper roll holder (rod) P-3050, a mobile drawer unit (similar to HB-580/3, page 113) and a swivel stool DH-2000 (page 104).

### SSL-6558

This examination unit includes a universal table (Article No. U-1065, page 49), two mountable melamine-coated base drawer units with two and three drawers respectively (similar to Article No. HB-580/2 and HB-580/3, page 113). The drawers units can be positioned individually as desired.

### UEH-5865

This unit differs from model SSL 6558 only in the nature of the drawer front finish. Here, the panels are covered with the same synthetic leather as the table.

Article Number	BUD-6265	SSL-6558	UEH-5865
Picture(s)	1	2	3
Price			
Length (mm) Headrest/Main section	1900 500/1400	1900 500/1400	
Width (mm)	650	650	
Height (mm)	650	650/700/750/800	
Frame	Lacquered	Lacquered	
Upholstery thickness (mm)	47	64	
Headrest adjustment	+30° (Grid)	+40° (Clamping rods)	
Maximum load (kg)	150	200	
Paper roll holder	✓	P-2065	
Accessories	—	Pages 60-61	

# Universal Tables

## Treatment Tables

The treatment tables RU-1062 and RUL-6000 stand out by their shapely design, enhanced by a tabletop slightly raised from the flat oval tubular steel leg and mattress frame by spacers, and finished with softly rounded corners. The treatment tables ULI and ELI are economical products which find their application in many areas.



1 RU-1062



2 RUL-6000



3 ULI-6000



4 ELI-8000

Article Number	RU-1062	RUL-6000	ULI-6000	ULI-6020	ULI-6030	ELI-8000	ELI-8020	ELI-8030
Picture(s)	1	2	3	Similar 3	Similar 3	4	Similar 4	Similar 4
Price								
Length (mm)	1920	1900		1900			1900	
Headrest/Main section	320/1600	500/1400		500/1400			500/1400	
Width (mm)	620	620		650			800	
Height (mm)	700	680		650			650	
Tabletop frame	Lacquered	Lacquered		Lacquered			Lacquered	
Table legs								
▪ Lacquered	✓	✓	✓	—	—	✓	—	—
▪ Chrome-plated	—	—	—	✓	—	—	✓	—
▪ Stainless steel	—	—	—	—	✓	—	—	✓
Upholstery thickness (mm)	64	47		47			47	
Headrest adjustment	-25° / +30° (Gas spring)	+30° (Grid)		+30° (Grid)			+30° (Grid)	
Maximum load (kg)	150	150		150			150	
Levelling feet	✓	✓		▪			▪	
Paper roll holder	PU-6000	PU-5500		PU-6000			PE-8000	
Accessories	—	—		—			Pages 60-61	

✓ included in delivery

▪ available (extra charge)

— not available



# Universal Tables

## Reclining Lounge Chair / EEG Chair

The reclining lounge chair, ideal for relaxing, also finds its use as EEG chair in the field of electroencephalography (EEG - a method for measuring electrical brain waves for medical diagnostic and neurological research). This model is particularly comfortable and user-friendly.



1 KL-1070/M, HRM-14/KL



2 KL-1070/M, HRM-14/KL

Article Number	KL-1070/M
Picture(s)	In 1-4
Price	
Length (mm) Back/Seat/Legs	1900 830/470/600
Width (mm)	580
Seat height (mm)	550
Dimensions when folded (mm)	1100x710x400
Tubular steel frame	Lacquered (Colour Basalt)
Upholstery (cover)	Synthetic leather
Upholstery thickness (mm)	40
Maximum load (kg)	140
Armrests	✓
Backrest/Leg section adjustment	Continuous
Accessories	Neck bolster with strap HRM-14/KL

### Additional Features

- This EEG chair is a 10 kg lightweight tubular steel construction.
- The leg section and the backrest can be adjusted together continuously via the durable plastic-coated footrest to achieve an ideal rest and relaxation position.
- The extra soft upholstery ensures maximum comfort during longer seating times. The padding is drawn over the covering and held by side spring clips.
- Wherever space is limited, the chair can be folded quickly and easily and stowed away when not in use.



3 KL-1070/M, HRM-14/KL



4 KL-1070/M - Detail

# Universal Tables

## Wall-Mounted Examination Tables

The wall-mounted examination tables are quick and easy to set up and intended for use in areas where a permanent table is not required or where space is limited. These tables are available in three frame finishes (lacquered, chrome, stainless steel) and two widths (650, 800 mm).



1 WKL-1065 - Headrest to the right



2 WKL-1065 - Headrest to the left

Article Number	WKL-1065	WKL-2065	WKL-3065	WKL-1080	WKL-2080	WKL-3080
Picture(s)	1, 2	Similar 1, 2	Similar 1, 2	Similar 1, 2	Similar 1, 2	Similar 1, 2
Price						
Frame	Lacquered	Chrome	Stainless steel	Lacquered	Chrome	Stainless steel
Length (mm) Headrest/Main section		1900 500/1400			1900 500/1400	
Length with wall brackets (mm)		2014			2014	
Width (mm)		650			800	
Height (mm)		700/750/800			900	
Upholstery thickness (mm)		64			64	
Headrest adjustment		+40° (Clamping rods)			+40° (Clamping rods)	
Maximum load (kg)		200			200	
Paper roll holder		PU-6000 (For wall mounting)			PE-8000 (For wall mounting)	
Accessories		Pages 60-61			Pages 60-61	

### Additional Features

- The desired position of the headrest, on the left or on the right-hand side, should be indicated at time of ordering.
- The individually foldable legs are equipped with a safeguard against accidental collapse.
- Reinforced wall brackets are included with the table.
- In the closed position, the distance from the tabletop to the wall is approximately 190 mm.

### Options

- In the standard version, these tables fold downwards, along the wall. This way, the visible part is the upholstered tabletop (not the leg frame with the underside of the table and the headrest adjustment mechanism).

The wall-mounted examination tables are also available:

- Folding upwards, horizontally, parallel to the wall (headrest left or right-hand side),
- Folding upwards, vertically, perpendicular to the wall (headrest towards the wall or towards the centre of the room), or
- With wall fixations stretching all the way down to the floor.

Please contact us for pricing information.

# Universal Tables

## Children / Pediatric Examination Tables

The pediatric examination tables have been especially designed for the treatment of children. The headrest and the leg section can be continuously adjusted up or down, or folded away completely to adapt to different children body sizes.

These children examination tables, with or without castors, are available in three frame finishes (lacquered, chrome-plated, stainless steel), two widths (650, 800 mm) as well as four heights (650, 700, 750, 800 mm).



1 KU-1065



2 KU-2065/LR

Article Number	KU-1065	KU-2065	KU-3065	KU-1080	KU-2080	KU-3080
Picture(s)	1	Similar 1	Similar 1	Similar 1	Similar 1	Similar 1
Price						
Frame	Lacquered	Chrome	Stainless steel	Lacquered	Chrome	Stainless steel
Length (mm) Head/Middle/Feet		1900 500/900/500			1900 500/900/500	
Width (mm)		650			800	
Height (mm)		650/700/750/800			650/700/750/800	
Tabletop		3 sections			3 sections	
Upholstery thickness (mm)		64			64	
Maximum load (kg)		80			80	
Head/Foot section adjustment		-90° / +30° (Clamping rods)			-90° / +30° (Clamping rods)	
Castors Ø 75 mm (two locking)	KU-1065/LR	KU-2065/LR	KU-3065/LR	KU-1080/LR	KU-2080/LR	KU-3080/LR
Paper roll holder	PU-6000 (To be mounted under the tabletop or on the wall)			PE-8000 (To be mounted under the tabletop or on the wall)		
Accessories	Pages 60-61			Pages 60-61		

## Additional Features

- The children examination tables can also be equipped with castors Ø 125 mm or with a central locking castor system. Pricing information is available upon request.
- Unless requested otherwise at time of ordering, the tables will be supplied in the standard height of 650 mm.

# Universal Tables

## Universal Tables with Continuous Headrest Adjustment

Universal tables can be described in many ways because of their versatility. Universal tables are also known as treatment couches, patient tables, examination tables, ECG tables, gymnastics tables, therapy tables, or as massage tables. The models presented here offer a wide range of headrest adjustment options as well as tabletop configurations.



1 U-1065




2 U-1065/3

Article Number	U-1065	U-2065	U-3065	EKG-1080	EKG-2080	EKG-3080
Picture(s)	1	3	Similar 3	Similar 1	Similar 3	Similar 3
Price						
Frame	Lacquered	Chrome	Stainless steel	Lacquered	Chrome	Stainless steel
Length (mm) Headrest/Main section		1900 500/1400			1900 500/1400	
Width (mm)		650			800	
Height (mm)		650/700/750/800			650/700/750/800	
Tabletop		2 sections			2 sections	
Upholstery thickness (mm)		64			64	
Maximum load (kg)		200			200	
Headrest adjustment		+40° (Clamping rods)			+40° (Clamping rods)	
Negative headrest adjustment +/-30°	SU-1065	SU-2065	SU-3065	SU-1080	SU-2080	SU-3080
Tabletop in 3 sections (500/500/900 mm) with foot section adjustment +20°	U-1065/3	U-2065/3	U-3065/3	EKG-1080/3	EKG-2080/3	EKG-3080/3
Paper roll holder		P-2065			P-2080	
Accessories		Pages 60-61			Pages 60-61	

### Additional Features

- The reinforced brackets in the frame corners and crossbars between the table legs ensure the stability and robustness of the structure.

- The standard height is 650 mm, unless requested otherwise upon ordering.
- The universal tables are also available with a width of 1000 mm (Pricing information upon request).

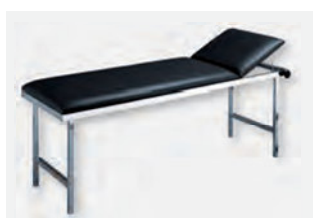

3 U-2065



4 SU-1065



# Universal Tables

## Universal Tables with Castors and Side Guards - Patient Transport Tables / Recovery Tables

The universal tables presented on the previous page are also available with castors as patient transport tables and with side guards for safe patient transport. They can therefore also be used as recovery tables.



1 U-1065/LR



2 AWL-1065, Korb

	Patient Transport Tables		Recovery Tables	
Article number	U-1065/LR	EKG-1080/LR	AWL-1065	AWL-1080
Picture(s)	1	Similar 1	In 2	Similar 2
Price				
Frame	Lacquered		Lacquered	
Length (mm)	1900		1900	
Headrest/Main section	500/1400		500/1400	
Width (mm)	650	800	650	800
Height (mm)	650/700/750/800		650/700/750/800	
Tabletop	2 sections		2 sections	
Upholstery thickness (mm)	64		64	
Maximum load (kg)	200		200	
Side guards (pivoting, 800x250 mm)	—		✓ (Overall width 795 mm)	✓ (Overall width 975 mm)
Castors	Ø 75 mm (two locking)		Ø 75 mm (two locking)	
Headrest adjustment	+40° (Clamping rods)		+40° (Clamping rods)	
Frame - Chrome	U-2065/LR	EKG-2080/LR	AWL-2065	AWL-2080
Paper roll holder	P-2065	P-2080	P-2065	P-2080
Storage basket	Korb		Korb	
Accessories	Pages 60-61		Pages 60-61	

✓ included in delivery — not available

### Additional Features

- These tables are also available with a width of 1000 mm. Further equipment include push bars (at the head and/or foot end), wall guards and an all-around frame protector. Please contact us for pricing information.
- The patient transport/recovery tables are also available with height adjustment to make getting on/off more convenient for the patient. These models are presented on pages 55-56.

- The tables described on this page are not suitable for magnetic resonance imaging (MRI). Please refer to pages 28-33 of this catalogue for MRI-compatible products.

### Accessories

- Wire storage basket (Article No. Korb), dimensions 550x290x190 mm (WxDxH), lacquered in light grey RAL 7035. Retrofitting is not possible.

# Universal Tables

## Patient Transport Table / Recovery Table with Central Locking Castors

This patient transport/recovery table stands out not only by its functionality but also by its user-friendly design. The two chrome-plated push bars and the large castors ensure excellent manoeuvrability and comfort of use during transport.

Article Number	PTL-1070/LRZ
Picture(s)	In 1-3
Price	
Length (mm)	1910 (Overall length 1950)
Headrest/Main section	560/1330
Width (mm)	700 (Overall width 750)
Height (mm)	650
Headrest adjustment	+30° (Gas pressure spring)
Upholstery thickness (mm)	64
Maximum load (kg)	200
Side guards - Pivoting	✓
Push bars	✓
Castors	Double castors Ø 125 mm with central locking system
Paper roll holder	P-3059/PTL
Storage basket	Korb
Infusion (IV) holder	INF-30170



1 PTL-1070/LRZ, INF-30170, Korb

### Additional Features

- All four central locking double castors are connected by rods. A foot pedal to secure/release is located at the head and foot ends of the table.
- All walls and doors as well as the patient transport table itself are protected by the black rubber guard mounted all around the tabletop frame and side guards.
- The headrest, main table surface and padded areas on the side guards are upholstered with synthetic leather.
- The side guards (760x230 mm, LxH) can be adjusted up or down with a simple rotation mechanism and locked in place with a sliding bolt. When both side guards are down, the overall width of the upholstered area is increased from 600 to 700 mm.
- The headrest is continuously adjustable up to +30° by gas pressure spring (300N). The release button is located at the front end, directly underneath the headrest (picture 3).

### Accessories

- At the front end (headrest), built-in brackets (Ø 18 mm), located on the right and left-hand sides are designed to hold the optional infusion (IV) holder. The infusion pole comes with a bottle holder and four hooks.
- The wire storage basket (550x290x190 mm, WxDxH), lacquered in light grey RAL 7035, can be pulled out from either side. Retrofitting is not possible.



2 PTL-1070/LRZ, INF-30170, Korb

### Options

- The patient transport and recovery tables are also available with height adjustment to make it more convenient for the patient to get on/off the table. These products are presented on pages 55-56.



3 PTL-1070/LRZ - Detail

# Universal Tables

## XXL Universal Tables for Bariatric Patients

We have extended our range of universal tables to include products particularly intended for bariatric patients. Adiposity (or obesity) is a nutritional and metabolic disease causing patient to be extremely overweight.

In addition to safety, these XXL tables provide patients with the necessary comfort during an examination or treatment.



1 XXL-1080/1900/PH

Article Number	XXL-1080/ 1900/PH	XXL-1080/ 2000/PH
Picture(s)	1-3	Similar 1-3
Price		
Frame	Lacquered	
Length (mm)	1900	2000
Headrest/Main section	500/1400	500/1500
Width (mm)	800	
Height (mm)	650	
Tabletop	2 sections	
Upholstery thickness (mm)	64	
Maximum load (kg)	300	
Headrest adjustment	+40° (Clamping rods)	
Paper roll holder (Rod at foot end)	✓	
Accessories	Pages 60-61	



2 XXL-1080/1900/PH - Detail



3 XXL-1080/1900/PH - Detail

### Options

- The XXL examination tables for bariatric patients are also available with hydraulic height adjustment (Article No. H-XXL 1080/Y) or electric height adjustment (Article No. H-XXL 1080/E). Please refer to the next page for more information.
- More XXL/bariatric products such as blood collection chairs are presented with the transfusion chairs on page 100.

### Accessories

- The XXL tables are also available with a face cutout (Article No. GS-1710) and upholstered insert (Article No. EfGS).
- Further accessories such as positioning aids, standard rails or side guards can be found on catalogue pages 60-61.

# Universal Tables

## XXL (Bariatric) Universal Tables with Hydraulic and Electric Height Adjustment

Our current XXL product range (examination tables with fixed height and blood collection chairs) now also includes two height-adjustable universal table versions (hydraulic and electric).

Both height-adjustable XXL tables, specifically intended for bariatric patients, can safely withstand a load up to 300 kg.

Article Number	H-XXL 1080/Y	H-XXL 1080/E
Picture(s)	1	2, 3
Price		
Frame	Lacquered	
Length (mm)	1900	
Headrest/Main section	500/1400	
Width (mm)	800	
Height adjustment (mm)	650-850 (Hydraulic)	650-850 (Electric)
Tabletop	2 sections	
Upholstery thickness (mm)	64	
Maximum load (kg)	300	
Headrest adjustment	+40° (Clamping rods)	
Levelling feet	✓	
Paper roll holder	P-2080/A	
Accessories	Pages 60-61	



1 H-XXL 1080/Y



2 H-XXL 1080/E



3 H-XXL 1080/E

### Additional Features

- The hydraulic height adjustment, without lateral offset, is done through a centrally built-in foot pump with little effort, even with heavier patients (up to 300 kg). The foot pedal can be mounted on either side (please specify upon ordering). A second foot pedal is available at extra charge (Article No. FPH-10).
- The electric height adjustment is done with a maintenance-free, spill-proof and extremely quiet electric motor, controlled by air pulses via a hand switch (or foot switch at no extra charge, please specify upon ordering). Drive specifications: 230V, 50/60Hz, 3,7A, 850W, IPX4, v=14mm/s, 7000N.
- The base of the frame is equipped with levelling feet to compensate for uneven floors and ensure a high stability.
- To prevent rocking and swinging movements of the tabletop, the XXL universal tables with hydraulic height adjustment are manufactured with precision dual-box rectangular profiles.

# Universal Tables

## Universal Tables with Manual Height Adjustment

To anyone looking for a universal table that only needs to be adjusted in height occasionally, we recommend a universal table with manual height adjustment as a cost-effective solution. The operation of the hand crank on this manually height-adjustable table is done from the foot end. After use, the crank handle can simply be pushed back under the tabletop.



1 H-U 1065/H



2 H-SU 1080/H

Article Number	H-U 1065/H	H-U 2065/H	H-U 3065/H	H-EKG 1080/H	H-EKG 2080/H	H-EKG 3080/H
Picture(s)	1	Similar 1	Similar 1	Similar 1	Similar 1	Similar 1
Price						
Frame	Lacquered	Chrome	Stainless steel	Lacquered	Chrome	Stainless steel
Length (mm) Headrest/Main section		1900 500/1400			1900 500/1400	
Width (mm)		650			800	
Height adjustment (mm)		650-850 (Manual)			650-850 (Manual)	
Tabletop		2 sections			2 sections	
Upholstery thickness (mm)		64			64	
Maximum load (kg)		200			200	
Headrest adjustment		+40° (Clamping rods)			+40° (Clamping rods)	
Negative headrest adjustment +/-30°	H-SU 1065/H	H-SU 2065/H	H-SU 3065/H	H-SU 1080/H	H-SU 2080/H	H-SU 3080/H
Paper roll holder		P-2065/A			P-2080/A	
Accessories		Pages 60-61			Pages 60-61	

## Additional Features

- The universal tables with manual height adjustment are manufactured with precision dual-box rectangular profiles to prevent rocking and swinging movements of the tabletop. The maintenance-free gear with optimal transmission ratio ( $i=10\text{mm/rev}$ ) ensures smooth, quick and easy operation of the hand crank.
- The chrome-plated crossbars fitted with levelling feet and the flat foot structure are more resistant than a lacquered frame against external influences such as repeated vacuum cleaner impacts.

- In the table, the term 'frame' only refers to the tabletop frame. The height adjustment mechanism is only available in lacquered finish and the leg frame (crossbars) in chrome finish.

## Options

- The paper roll holder for the H-SU... models may be ordered under Article No. PU-6000 (for table width 650 mm) and PE-8000 (for table width 800 mm).
- The universal tables with manual height adjustment are also available with a width of 1000 mm, or with a three-piece tabletop. Pricing information is available upon request.



# Universal Tables

## Universal Tables with Hydraulic Height Adjustment

Another type of height adjustment for the universal tables is a hydraulic system, operated from a centrally built-in foot pump. The height adjustment can be performed with little effort, even with heavier patients (up to 200 kg). The foot pedal can be mounted on either side (please specify upon ordering). A second foot pedal is available at extra charge.



1 H-U 1065/Y



2 H-SU 1080/Y

Article Number	H-U 1065/Y	H-U 2065/Y	H-U 3065/Y	H-EKG 1080/Y	H-EKG 2080/Y	H-EKG 3080/Y
Picture(s)	1	Similar 1	Similar 1	Similar 1	Similar 1	Similar 1
Price						
Frame	Lacquered	Chrome	Stainless steel	Lacquered	Chrome	Stainless steel
Length (mm) Headrest/Main section		1900 500/1400			1900 500/1400	
Width (mm)		650			800	
Height adjustment (mm)		650-850 (Hydraulic)			650-850 (Hydraulic)	
Tabletop		2 sections			2 sections	
Upholstery thickness (mm)		64			64	
Maximum load (kg)		200			200	
Headrest adjustment		+40° (Clamping rods)			+40° (Clamping rods)	
Negative headrest adjustment +/-30°	H-SU 1065/Y	H-SU 2065/Y	H-SU 3065/Y	H-SU 1080/Y	H-SU 2080/Y	H-SU 3080/Y
Add. foot pedal		FPH-10			FPH-10	
Paper roll holder		P-2065/A			P-2080/A	
Accessories		Pages 60-61			Pages 60-61	

### Additional Features

- The precision dual-box rectangular profiles help prevent rocking and swinging movements of the tabletop.
- The chrome-plated crossbars fitted with levelling feet are more durable against external influences such as repeated vacuum cleaner impacts.
- In the table, the term 'frame' only refers to the tabletop frame. The height adjustment mechanism is only available in lacquered finish and the leg frame (crossbars) in chrome finish.

### Options

- The paper roll holder for the H-SU... models may be ordered under Article No. PU-6000 (for table width 650 mm) and PE-8000 (for table width 800 mm).
- The universal tables with hydraulic height adjustment are also available with a width of 1000 mm or with a three-piece tabletop. Please contact us for pricing information.

# Universal Tables

## Universal Tables and Recovery Tables with Hydraulic Height Adjustment

The height-adjustable universal tables with scissor mechanism are characterized by their particularly large adjustment range. The hydraulic height adjustment (500-980 mm) is done by a foot pump. For added convenience, the positioning of the foot pedal is possible on the left or the right-hand side.



1 H-U 1065/Y 980



2 H-AWL 1065/Y 980

Article Number	H-U 1065/Y 980	H-EKG 1080/Y 980	H-AWL 1065/Y 980	H-AWL 1080/Y 980
Picture(s)	1	Similar 1	2	Similar 2
Price				
Frame	Lacquered		Lacquered	
Length (mm)	1900		1900	
Headrest/Main section	500/1400		500/1400	
Width (mm)	650	800	650	800
Height adjustment (mm)	500-980 (Hydraulic)		500-980 (Hydraulic)	
Tabletop	2 sections		2 sections	
Upholstery thickness (mm)	64		64	
Maximum load (kg)	200		200	
Headrest adjustment	+40° (Clamping rods)		+40° (Clamping rods)	
Central locking system	▪		✓ (Castors Ø 125 mm)	
Side guards (pivoting)	▪		✓ (Overall width 795 mm)	✓ (Overall width 975 mm)
Frame - Chrome	H-U 2065/Y 980	H-EKG 2080/Y 980	H-AWL 2065/Y 980	H-AWL 2080/Y 980
Add. foot pedal	FPH-10		FPH-10	
Shock position/ Trendelenburg 11°	SL-AWL		SL-AWL	
Paper roll holder	P-2065/A	P-2080/A	P-2065/A	P-2080/A
Accessories	Pages 60-61		Pages 60-61	

✓ included in delivery    ▪ available (extra charge)

### Additional Features

- In the table, 'frame' only refers to the tabletop frame (with side guards, 760x250 mm, LxH). The height adjustment mechanism and leg frame are available in lacquered finish only.
- The products described on this page are not suitable for magnetic resonance imaging (MRI). The MRI patient transport tables can be found on pages 28-33 of this catalogue.

### Options

- The red release button for the shock position is embedded in the tabletop frame at the foot end. Retrofitting is not possible.
- The patient transport/recovery tables can be equipped with push bars, wall guards and all-around frame protector. Please contact us for pricing information.

# Universal Tables

## Recovery Tables with Electric and Battery Powered Height Adjustment

We have expanded our product range of recovery/patient transport tables to electric height adjustable and battery powered versions. These tables are equipped with a central locking castor system and two side guards.



1 H-AWL 1080/E 930



2 H-AWL 1065/AC 1070

Article Number	H-AWL 1065/E 930	H-AWL 1080/E 930	H-AWL 1065/AC 1070	H-AWL 1080/AC 1070
Picture(s)	Similar 1	1	2	Similar 2
Price				
Frame	Lacquered		Lacquered	
Length (mm)	1900		1900	
Headrest/Main section	500/1400		500/1400	
Width (mm)	650	800	650	800
Height adjustment (mm)	500-930 (Electric)		670-1070 (Battery)	
Tabletop	2 sections		2 sections	
Upholstery thickness (mm)	64		64	
Maximum load (kg)	200		150	
Headrest adjustment	+40° (Clamping rods)		+40° (Clamping rods)	
Central locking system	✓ (Castors Ø 125 mm)		✓ (Castors Ø 125 mm)	
Side guards (pivoting)	✓ (Overall width 795 mm)	✓ (Overall width 975 mm)	✓ (Overall width 795 mm)	✓ (Overall width 975 mm)
Frame - Chrome	H-AWL 2065/E 930	H-AWL 2080/E 930	H-AWL 2065/AC 1070	H-AWL 2080/AC 1070
Shock position/ Trendelenburg 11°	SL-AWL		SL-AWL	
Battery module	-		BMEAC	
Paper roll holder	P-2065/A	P-2080/A	P-2065/A	P-2080/A
Accessories	Pages 60-61		Pages 60-61	

### Additional Features

- In the table, 'frame' only refers to the tabletop frame (with side guards). The height adjustment mechanism and leg frame are available in lacquered finish only.
- Drive specifications (motor):  
230V, 50/60Hz, 3,7A, 850W, IPX4, v=14mm/s, 7000N.  
Drive specifications (each battery unit):  
230V, 50/60Hz, 4,0A, 160W, IPX6, v=16mm/s, 2x2000N.

### Options

- The red release button for the shock position on the electric adjustable patient transport tables is embedded in the tabletop frame at the foot end. Retrofitting is not possible.
- Push bars, wall guards and all-around frame protector are available. Please contact us for pricing information.

# Universal Tables

## Universal Tables with Electric Height Adjustment

Besides the manual and hydraulic height-adjustable universal tables, we also offer models with electric height adjustment. The universal tables with scissor mechanism are characterized by their particularly large height adjustment range of 430 mm.

The models shown here are available in a variety of frame finishes, widths and adjustment options for the headrest.



1 H-U 1065/E



2 H-U 1065/E 930

Article Number	H-U 1065/E	H-EKG 1080/E	H-U 1065/E 930	H-EKG 1080/E 930
Picture(s)	1	Similar 1	2	Similar 2
Price				
Frame	Lacquered		Lacquered	
Length (mm) Headrest/Main section	1900 500/1400		1900 500/1400	
Width (mm)	650	800	650	800
Height adjustment (mm)	650-850 (Electric)		500-930 (Electric)	
Tabletop	2 sections		2 sections	
Upholstery thickness (mm)	64		64	
Maximum load (kg)	200		200	
Headrest adjustment	+40° (Clamping rods)		+40° (Clamping rods)	
Frame - Chrome	H-U 2065/E	H-EKG 2080/E	H-U 2065/E 930	H-EKG 2080/E 930
Frame - Stainless steel	H-U 3065/E	H-EKG 3080/E	—	—
Negative headrest adjustment +/-30°	H-SU 1065/E	H-SU 1080/E	—	—
Paper roll holder	P-2065/A	P-2080/A	P-2065/A	P-2080/A
Accessories	Pages 60-61		Pages 60-61	

### Additional Features

- The height adjustment is done with a maintenance-free, spill-proof and extremely quiet electric motor, controlled by air pulses via a hand switch (or foot switch at no extra charge, please specify upon ordering). Drive specifications: 230V, 50/60Hz, 3,7A, 850W, IPX4, v=14mm/s, 7000N.
- In the table, 'frame' refers to the tabletop frame only.

### Options

- The universal tables with electric height adjustment are also available with a width of 1000 mm or with a 3-section tabletop. Pricing information is available upon request.



3 H-SU 1065/E - Detail

# Universal Tables

## AGA-DUO-LIFT Universal Tables with Electric Height and Tilt Adjustment

A further version of the electric height-adjustable universal tables is the AGA-DUO-LIFT. These tables offer, in the standard version, electric tilt adjustment in addition to electric height adjustment.

The AGA-DUO-LIFT universal tables are available in three frame finishes (lacquered, chrome, stainless steel) and two widths (650, 800 mm).



1 H-U 1065/EE 850



2 H-U 1065/EE 850

Article Number	H-U 1065/ EE 850	H-U 2065/ EE 850	H-U 3065/ EE 850	H-EKG 1080/ EE 850	H-EKG 2080/ EE 850	H-EKG 3080/ EE 850
Picture(s)	1, 2	Similar 1, 2	Similar 1, 2	Similar 1, 2	Similar 1, 2	Similar 1, 2
Price						
Frame	Lacquered	Chrome	Stainless steel	Lacquered	Chrome	Stainless steel
Length (mm) Headrest/Main section	1900 500/1400			1900 500/1400		
Width (mm)	650			800		
Height adjustment (mm)	600-850			600-850		
Adjustment mechanism	Electric, with hand switch (6 buttons for simultaneous height and tilt adjustment)			Electric, with hand switch (6 buttons for simultaneous height and tilt adjustment)		
Trendelenburg/ Reverse Trendelenburg	11°			11°		
Tabletop	2 sections			2 sections		
Upholstery thickness (mm)	64			64		
Maximum load (kg)	150			150		
Headrest adjustment	+40° (Clamping rods)			+40° (Clamping rods)		
Paper roll holder	P-2065/A			P-2080/A		
Accessories	Pages 60-61			Pages 60-61		

### Additional Features

- In the table the term 'frame' only refers to the tabletop frame. The leg frame is only available in lacquered finish.
- The smooth operating anodized light alloy lifting columns (low-voltage electric drives) are controlled via a transformer by a hand module with six buttons allowing simultaneous height and tilt adjustment (Multitasking).  
Drive specifications (each motor unit): 230V, 50/60Hz, 6,0A, 240W, IPX6, v=16mm/s, 2x1000N.
- Please note that in a slightly modified form, the universal tables described on this page are also available as X-ray tables, for use with a C-arm. Please see pages 22-25.

### Options

- The universal tables with electric height adjustment are also available with a 3-section tabletop. Pricing information is available upon request.










# Accessories

## Accessories for Universal Tables

Picture	Article Number	Description	Price
	RLM-20	Knee roll, Ø 180x560 mm, synthetic leather, with zipper closure. Upholstery colours on page 4.	
	RM-12	Neck roll, Ø 140x400 mm, synthetic leather, with zipper closure. Upholstery colours on page 4.	
	KKM-16	Wedge-shaped headrest, dimensions 450x400x80/10 mm (LxWxH), synthetic leather, with zipper closure. Upholstery colours on page 4.	
	HRM-14	Half-round neck bolster, Ø 200x400 mm, synthetic leather, with zipper closure. Upholstery colours on page 4.	
	HRM-14/Band HRM-14/Klett	Half-round neck bolster with strap. Half-round neck bolster with hook-and-loop fasteners.  Synthetic leather, Ø 200x400 mm. Upholstery colours on page 4.	
	PU-6000 PE-8000 PI-10000	Paper roll holder, chrome-plated, width 650 mm. Paper roll holder, chrome-plated, width 800 mm. Paper roll holder, chrome-plated, width 1000 mm.  Suitable for wall-mounting and for AGA treatment tables as well as AGA tables with negative headrest adjustment (SU models).	
	P-2065 P-2080 P-20100	Paper roll holder, chrome-plated, width 650 mm. Paper roll holder, chrome-plated, width 800 mm. Paper roll holder, chrome-plated, width 1000 mm.  Suitable for mounting on the leg frame of AGA tables with fixed height.	
	P-2065/A P-2080/A P-20100/A	Paper roll holder, chrome-plated, width 650 mm. Paper roll holder, chrome-plated, width 800 mm. Paper roll holder, chrome-plated, width 1000 mm.  With mounting hardware, to be attached on the tabletop frame of the height-adjustable AGA tables (except SU models).	
	LR-75 LR-125	Set of castors, Ø 75 mm. Set of castors, Ø 125 mm.  A set consists of four castors, two of which with locking function. At time of ordering, please indicate the required total height of the table (including castors).	
	LR/F-75 LR/F-125	Mobile base on retractable castors, Ø 75 mm. Mobile base on retractable castors, Ø 125 mm.  The foot pedal to secure/release is located at the head and foot end (press twice). At time of ordering, please indicate the required total height of the table (including castors).  *Retrofitting is not possible.	

## Accessories for Universal Tables

Picture	Article Number	Description	Price
	LR/Z-125	<p>Central locking castor system, Ø 125 mm, offering excellent mobility and comfort of use. All four double castors are connected together by rods. A foot pedal to secure/release is located at the head and foot ends of the table (press once).</p> <p>At time of ordering, please indicate the required total height of the table (including castors).</p> <p>*Retrofitting is not possible.</p>	
	SL-AWL	<p>Gas pressure spring at the foot end of the tabletop frame for the shock/Trendelenburg position (11°). The red release button is embedded in the tabletop frame.</p> <p>*Retrofitting is not possible.</p>	
	KVG	<p>Headrest adjustment by gas pressure spring instead of grid or clamping rod mechanism, for most AGA universal tables.</p> <p>*Retrofitting is not possible.</p>	
	GS-1710	<p>Face cutout in the headrest, oval shape, dimensions 205x105 mm (LxW). *Retrofitting is not possible.</p>	
	EfGS	<p>Upholstered insert to close the face cutout.</p>	
	PG-220 PG-260	<p>Patient strapping belt, dimensions 2200x250 mm (LxW). Patient strapping belt, dimensions 2600x250 mm (LxW).</p> <p>To secure the patient's position. Upholstery colours on page 4.</p>	
	OPS-40/A OPS-60/A OPS-100/A OPS-140/A OPS-170/A	<p>Standard rails (25x10 mm, HxW), set of two, stainless steel, with spacers. Available in different lengths (the digits in the article number indicate the respective length in cm).</p> <p>Please specify the desired position of the rail on the table frame (head, foot or center of the table) at time of order.</p> <p>*Retrofitting is possible to a limited extent only.</p>	
	SG-1000 SG-2000 SG-3000	<p>Side guard (removable), lacquered, dimensions 800x250 mm (LxH). Side guard (removable), dimensions 800x250 mm (LxH). Side guard (removable), dimensions 800x250 mm (LxH).</p> <p>*Retrofitting is not possible.</p>	
	SG-1000/1 SG-2000/1 SG-3000/1	<p>Side guard (pivoting), lacquered, dimensions 800x250 mm (LxH). Side guard (pivoting), chrome-plated, dimensions 800x250 mm (LxH). Side guard (pivoting), stainless steel, dimensions 800x250 mm (LxH).</p> <p>*Retrofitting is possible to a limited extent only.</p>	
	INF-1720/VA	<p>Infusion (IV) holder, stainless steel, length 1700 mm, with four hooks and fixing clamp (screw system). The clamp is fixed in the frame between the headrest and the main tabletop.</p> <p>*Retrofitting is possible to a limited extent only.</p>	